

AD-A129 961

SEPARATION OF ORGANOPHOSPHATE COMPOUNDS BY DYE-ASSISTED
CHROMATOGRAPHY(U) ARIZONA UNIV TUCSON DEPT OF CHEMISTRY
T GNANASAMBANDAN ET AL. JUN 83 N00014-81-K-0576


1/1

UNCLASSIFIED

F/G 7/4

NL



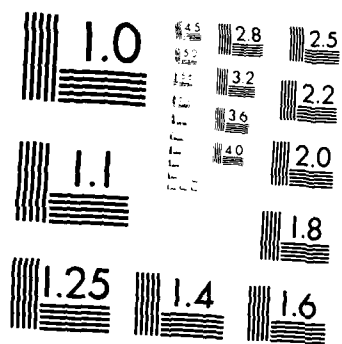

MICROCOPY RESOLUTION TEST CHART
NATIONAL BUREAU OF STANDARDS 1963-A

REPORT DOCUMENTATION PAGE

READ INSTRUCTIONS
BEFORE COMPLETING FORM

1. REPORT NUMBER 7	2. GOVT ACCESSION NO. ADA129961	3. RECIPIENT'S CATALOG NUMBER
4. TITLE (and Subtitle) "Separation of Organophosphate Compounds by Dye-Assisted Chromatography"		5. TYPE OF REPORT & PERIOD COVERED Technical & Interem
AUTHOR(s) T. Gnanasambandan and H. Freiser		6. PERFORMING ORG. REPORT NUMBER
PERFORMING ORGANIZATION NAME AND ADDRESS Department of Chemistry University of Arizona Tucson, Arizona 85721		8. CONTRACT OR GRANT NUMBER(s) N00014-81-K-0576
CONTROLLING OFFICE NAME AND ADDRESS Materials Sciences Division Office of Naval Research Arlington, VA 22217		10. PROGRAM ELEMENT, PROJECT, TASK AREA & WORK UNIT NUMBERS
4. MONITORING AGENCY NAME & ADDRESS (if different from Controlling Office)		12. REPORT DATE June, 1983
		13. NUMBER OF PAGES 11
		15. SECURITY CLASS. (of this report) unclassified
		15a. DECLASSIFICATION/DOWNGRADING SCHEDULE

16. DISTRIBUTION STATEMENT (of this Report)
This document has been approved for public release and sale; its distribution is unlimited.

17. DISTRIBUTION STATEMENT (of the abstract entered in Block 20, if different from Report)

18. SUPPLEMENTARY NOTES

19. KEY WORDS (Continue on reverse side if necessary and identify by block number)

dye, reversed, liquid chromatography, organophosphates, Brilliant Green, phosphonates, phosphites, phosphates

20. ABSTRACT (Continue on reverse side if necessary and identify by block number)

Dye-assisted, reversed phase, liquid chromatography has been applied to the separation of trace levels of a series of organophosphate compounds. Using dioxane-water mixtures containing $10^{-4}M$ Brilliant Green as mobile phases and either Partisil-ODS or Lichrosorb-ODS columns. The chromatographic separation of phosphonates, phosphites, and phosphates were accomplished.

DTIC
ELECTE
JUN 29 1983
S E D

ADA129961

DTIC FILE COPY

OFFICE OF NAVAL RESEARCH
Contract N00014-81-K-0576

TECHNICAL REPORT NO. 7

"Separation of Organophosphate Compounds
by Dye-Assisted Chromatography"

by

T. Gnanasambandan and H. Freiser

University of Arizona
Department of Chemistry
Tucson, Arizona 85721

Accession For	
NTIS GRA&I	<input checked="checked" type="checkbox"/>
DTIC TAB	<input type="checkbox"/>
Unannounced	<input type="checkbox"/>
Justification	
By	
Distribution/	
Availability Codes	
Dist	Avail and/or Special
A	

June, 1983

Reproduction in whole or in part is permitted for
any purpose of the United States Government

This document has been approved for public release
and sale; its distribution is unlimited

Separation of Organophosphate Compounds

Organophosphates represent a group of compounds that have a wide variety of applications ranging from chemical agents to pesticides that require a simple selective and sensitive trace level analytical technique.

Earlier efforts to develop an analytical method for these compounds utilized the more common chromatographic techniques. Thin layer and liquid chromatography were found to be inadequate because of their lack of sensitivity.(1) GLC with flame photometric detection of phosphorus has been used in the analyses of these compounds but some of these compounds need prior derivatization and most of the compounds cannot be separated with a single column. Ion chromatography(2) offers a sensitive technique for compounds that are phosphoric acids. The esters could be hydrolyzed to the corresponding phosphoric acid for determination. Compounds of similar pKa are found to be difficult to separate, however.

In our last publication(3) we reported a novel HPLC analytical technique for detecting and separating neutrals that are non-chromophoric. The purpose of this report is to show the use of the above technique in the separation of organophosphate compounds and also to establish the chromatographic conditions, that can be used for trace analysis.

Experimental

Materials. The organophosphates used for the study were obtained from various sources and were sometimes purified before

use. Degassed and deionized water, dioxane and methanol were used for all mobile phase preparations. Brilliant Green (Baker Analyzed) was used as received.

Purification of organophosphates. A 1.0ml sample is diluted to 10mL with Analar grade hexane. This solution is twice extracted with 10mL of 1% K_2CO_3 . The organic layer is separated and dried over Na_2SO_4 (anhydrous) and 1mL aliquots of the organic phase are diluted with methanol to 10 mL MeOH for analysis. It is necessary to prepare fresh samples to the dilution desired, because the solutions do deteriorate within two days.

Apparatus: A modular, high performance liquid chromatograph consisting of an Altex pump (Model 110A), a variable wavelength UV visible detector (Schoeffel Model SF770 Spectraflow), a refractive index detector (Showa Denko KK Shodex RI model SE-11), and a strip chart recorder (Linear Instruments model 261 m/m) was used. Sample injector was 10 uL loop injector (Spectra Physics rotary valve injector).

Mobile phase. The mobile phases used in this study were either dioxane/water or methanol/water. Each mobile phase also contained $1 \times 10^{-4} M$ Brilliant Green as the third component.

Columns. Two columns were used in the study. (a) Partisil ODS-2 10u 25cmx0.04cm and (b) Lichrosorb ODS-10u 25cmx0.04cm. Both these columns were slurry packed in our laboratory and had measured efficiencies of about 16,000 plates/meter. The cleaning and determination of the efficiencies of the column are determined as described elsewhere.(3)

Results and Discussion

Several investigators(4) have explored various chromatographic techniques for separating the major organophosphate esters. The present study is different in that we used a dye (Brilliant Green) as a component of the mobile phase to form an 'associate' with the analytes to enhance detection and separation.

Table 1 shows the various organophosphates chosen for our study. Ten compounds were chosen as a representative of the four families--phosphonates, phosphites, phosphate and phosphoramidate with a molecular weight range from 113-265.

The first set of experiments were conducted with a Lichrosorb ODS column using a mobile phase containing a mixture of dioxane and water. Dioxane/water mixture was chosen over methanol/water because our earlier studies indicated more efficient separation could be achieved. As we increased the organic modifier, however, dioxane had to be replaced with methanol to avoid creating a high back pressure resulting from the higher viscosity of dioxane.

Table 2 shows the capacity factors obtained for these compounds with various mobile phases. The data shows that good resolution was obtained for these compounds. Using a refractive index detector, we also obtained the capacity factors for these compounds on a dye-free column and using a dye-free mobile phase.

As in our earlier study with alcohols, the use of the dye increased analyte retention, a further indication of dye-analyte interaction.

Figure 1 shows the calibration curve obtained for these compounds. Peak height was plotted as a function of weight of organophosphate. As indicated from the graph, good linearity was obtained in the microgram range.

In separating, tributyl phosphate, dibutyl phosphite and dibutyl butyl phosphonate we had to use 70% methanol/water containing 10^{-4} M Brilliant Green as the mobile phase. Table 3 lists the various capacity factors obtained. As observed earlier and seen here as well, an increase of organic modifiers in the mobile phase will also result in the decrease of dye held on the column. Usually a negative peak with a characteristic capacity factor was seen after the samples were eluted. Table 3 also indicates that a separation between tributyl phosphate and dibutyl butyl phosphonate is not feasible. Figure 3 shows the calibration curves obtained for these compounds. They show good linearity and could easily be used for analytical purposes.

Although we did not use an extensive number of compounds to study in detail the mechanism of separation we have here clearly established that separation between the organophosphates is feasible and at an order of magnitude better sensitivity than obtained with normal modes of chromatography. Further work is under way.

References

1. Crabtree, E. V., and E. W. Sarver, "Review of Analytical Procedures for GB, VX, and their Degradation Products." Edgewood Arsenal Special publication EC-SP-76021, January 1977.
2. Schiff et al., "Ion Chromatographic Analysis of Environmental Pollutants," J.D. Mulik & E. Sawicki, Eds., Ann Arbor Sc., 2, 329,
3. H. Freiser and T. Gnanasambandan, Anal. Chem., 54, 1282-1285, 1982.
4. Akoi, Y., M. Takeda, and M. Uchiyama, J. Assoc. Offic. Anal. Chem. 58, 1286, 1975.

Legend for Figures

Figure 1. Separation of Organophosphates.

- 1.) Diethylphosphite
- 2.) Hexamethylphosphoramide
- 3.) Diethylethylphosphonate

Sample size 10uL contains about 80ug.

Experimental conditions: Column Spherisorb ODS 5u 25cm x 4.0mm

Flowrate: 1 mL/min. Detector wavelength: 436nm.

Figure 2. Calibration waves of various organophosphates.

Experimental conditions: Column Lichrosorb 10u ODS 25cm x 4.0mm

Mobile phase Dioxane/water 10:90 V/V% with 10^{-4} M Brilliant Green. Flow rate 1 mL/min. detection wavelength 436nm.

Figure 3. Calibration curve of various organophosphosphate.

Experimental conditions: Column Partisil ODS-2-10uM (25cmx0.04)

Mobile phase 75:25 V/V% MeOH/H₂O with 10^{-4} M Brilliant Green. Flow rate 1 mL/min. detection at 620nm

Table 1

Capacity Factors of Organophosphates (I)

MOBILE PHASE	30% MeOH/H ₂ O V/V WITH 10 ⁻⁴ M BG		10% DIOXANE/H ₂ O WITH 10 ⁻⁴ M BG
	ANALYTE	k'	k'
(DMMP)	Dimethylmethyl phosphonate	.6	.8
(DEEP)	Diethylethyl Phosphonate	4.4	6.9
(HMNP)	Hexamethyl Phosphoramidate	3.1	3.5
(DEP)	Diethyl Phosphite	1.6	2.5
	Phosdrin (methyl 3-hydroxy and crotonate, dimethyl phos- phate)	6.3	ND
(DDVP)	Dichlorvos	2.8	ND
	(2,2,-dichlorovinyl dimethyl phosphate)		

Column Partisil ODS-2 10um (25cm x 4.0mm id)

Flowrate 1nl per min.

Detector visible 436

Table 2

Capacity of Factors of Organophosphates (II)

MOBILE PHASE 75% MeOH/H₂O V/V WITH 10⁻⁴ B.GREEN

ANALYTE		k'
(DBP)	Dibutyl Phosphite	1.3
(DBBP)	Dibutylbutyl Phosphonate	3.2
(TBP)	Tributyl Phosphate	3.2
(TTP)	Tritolyl Phosphate	5.6

Column Partisil ODS-2 10u (25cm x 4.0mm id)

Flowrate 1 mL/min.

Detector wavelength 625nm.

Figure 1

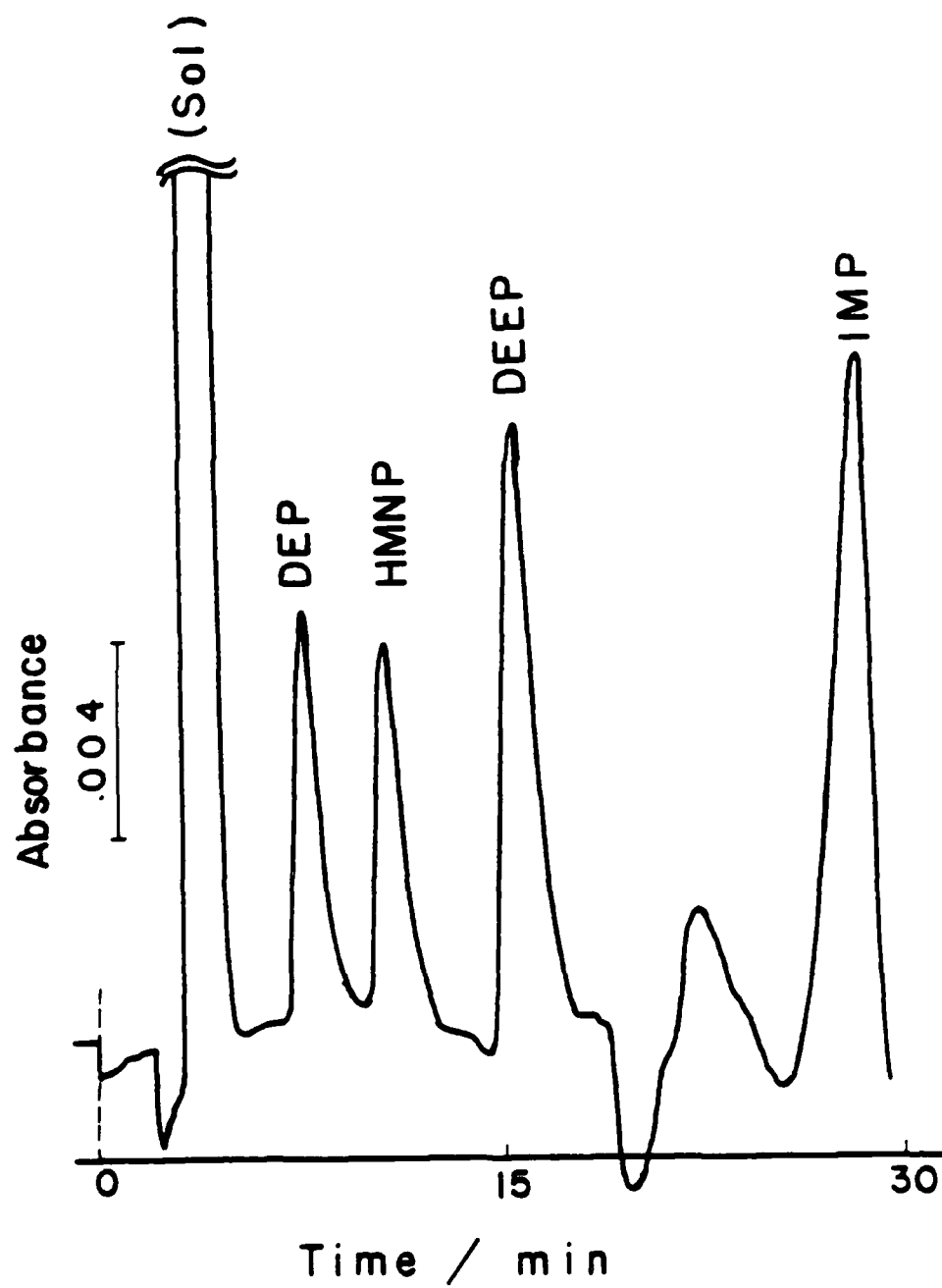


Figure 2


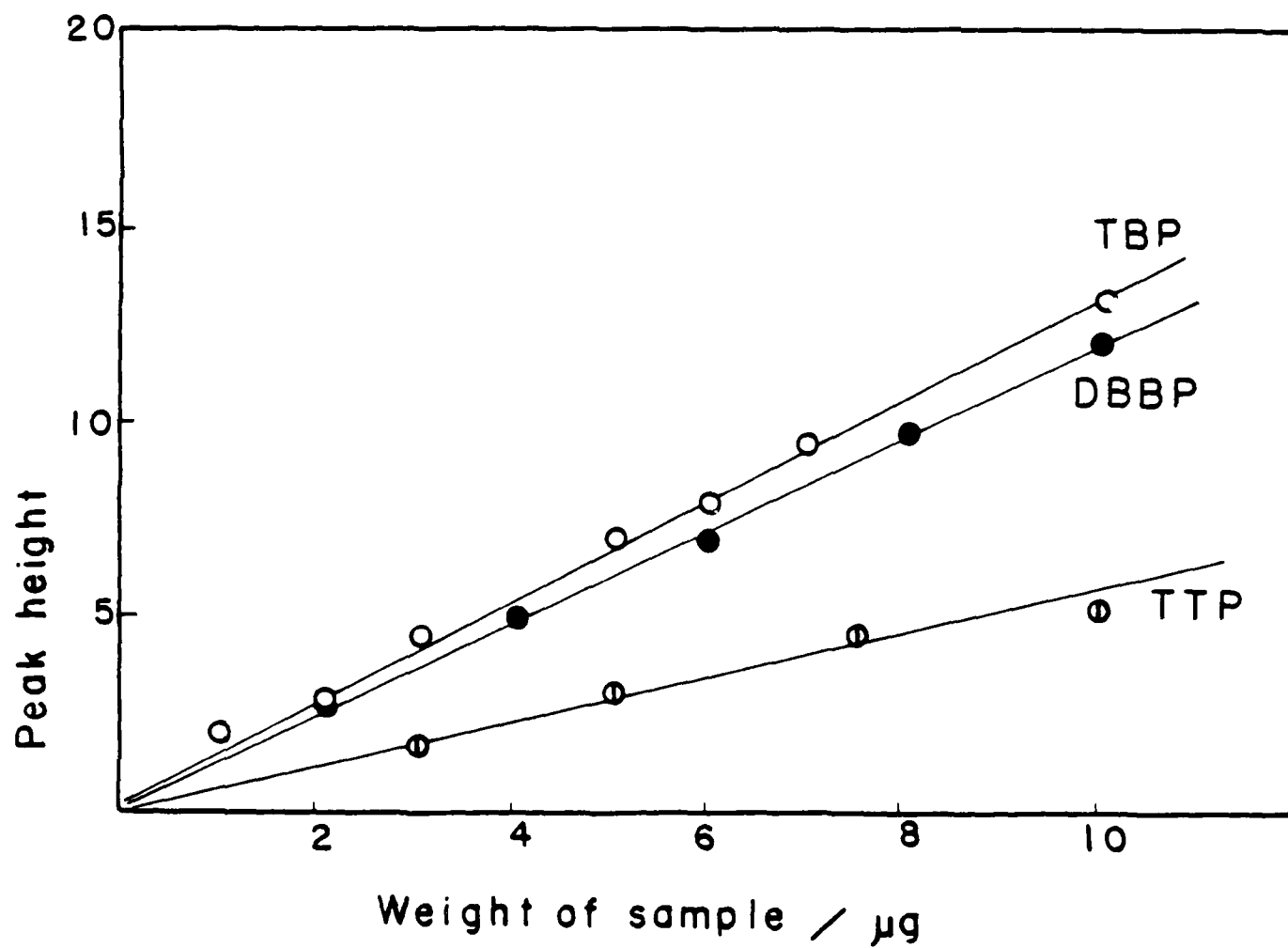


Figure 3

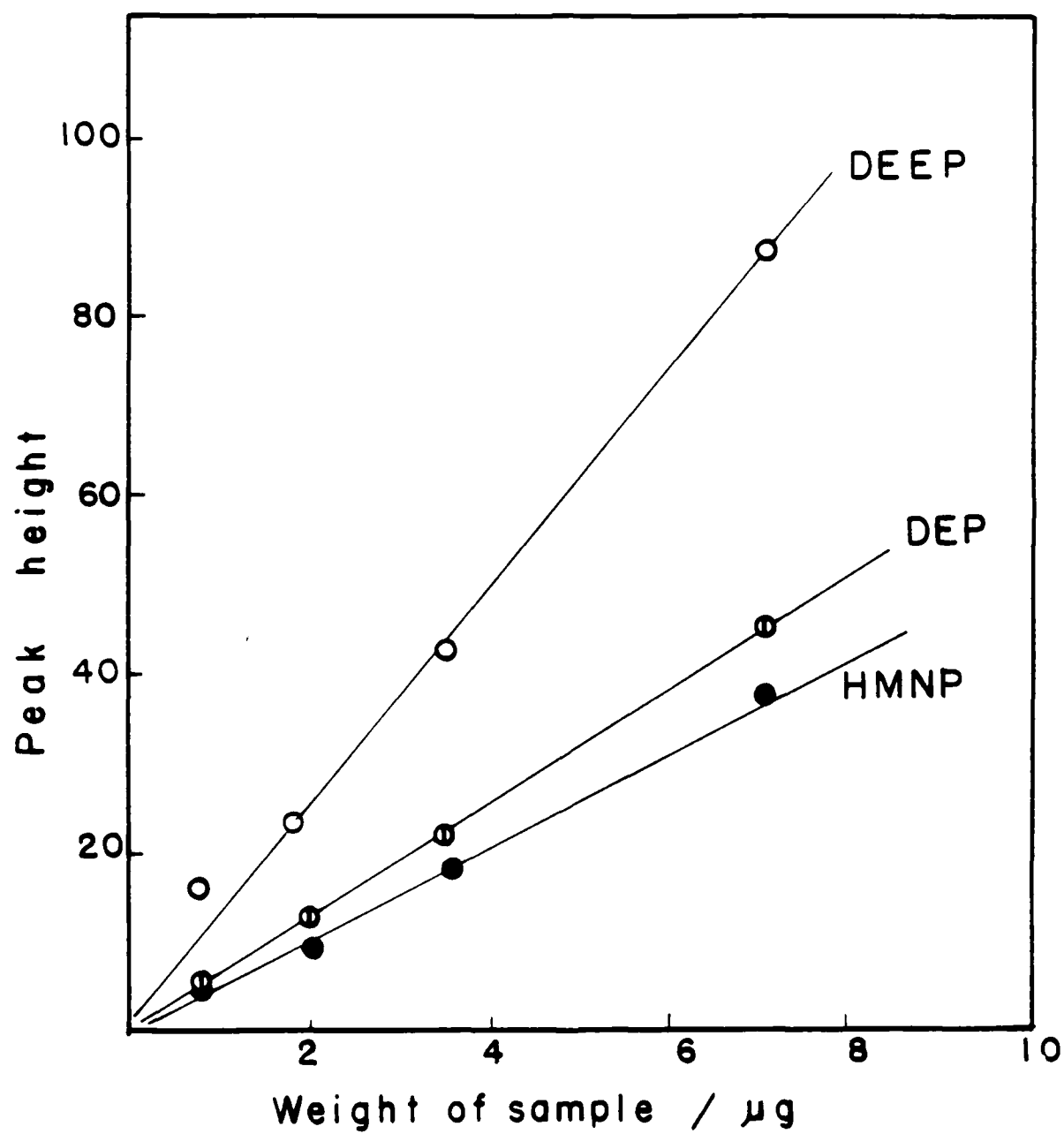

TECHNICAL REPORT DISTRIBUTION LIST, 051C

N
Cop

	<u>No.</u> <u>Copies</u>	
Dr. M. B. Denton Department of Chemistry University of Arizona Tucson, Arizona 85721	1	Dr. L. Jarris Code 6100 Naval Research Laboratory Washington, D.C. 20375
Dr. R. A. Osteryoung Department of Chemistry State University of New York at Buffalo Buffalo, New York 14214	1	Dr. John Duffin, Code 62 Dn United States Naval Postgraduate School Monterey, California 93940
Dr. J. Osteryoung Department of Chemistry State University of New York Buffalo, New York 14214	1	Dr. G. M. Hieftje Department of Chemistry Indiana University Bloomington, IN 47401
Dr. B. R. Kowalski Department of Chemistry University of Washington Seattle, Washington 98105	1	Dr. Victor L. Rehn Naval Weapons Center Code 3813 China Lake, California 93555
	1	Dr. Christie G. Enke Michigan State University Department of Chemistry East Lansing, Michigan 48824
Dr. D. L. Venezky Naval Research Laboratory Code 6130 Washington, D.C. 20375	1	Dr. Kent Eisentraut, MBT Air Force Materials Laboratory Wright-Patterson AFB, Ohio 45433
Dr. H. Freiser Department of Chemistry University of Arizona Tucson, Arizona 85721		Walter G. Cox, Code 3632 Naval Underwater Systems Center Building 148 Newport, Rhode Island 02840
Dr. H. Chernoff Department of Mathematics Massachusetts Institute of Technology Cambridge, Massachusetts 02139	1	Professor Isiah M. Warner Department of Chemistry Emory University Atlanta, Georgia 30322
Dr. A. Zirino Naval Undersea Center San Diego, California 92132	1	Professor George H. Morrison Department of Chemistry Cornell University Ithaca, New York 14853

TECHNICAL REPORT DISTRIBUTION LIST, 051C

	<u>No.</u> <u>Copies</u>
Professor J. Janata Department of Bioengineering University of Utah Salt Lake City, Utah 84112	1
Dr. Carl Heller Naval Weapons Center China Lake, California 93555	1
Dr. Denton Elliott AFOSR/NC Bolling AFB Washington, D.C. 20362	
Dr. J. Decorpo NAVSEA-05R14 Washington, D.C. 20362	
Dr. B. E. Spielvogel Inorganic and Analytical Branch P. O. Box 12211 Research Triangle Park, NC 27709	
Dr. Charles Anderson Analytical Chemistry Division Athens Environmental Lab. College Station Road Athens, Georgia 30613	
Dr. Samuel P. Perone L-326 LLNL Box 808 Livermore, California 94550	
Dr. B. E. Douda Chemical Sciences Branch Code 4052 Naval Weapons Support Center Crane, Indiana 47522	
Ms. Ann De Witt Material Science Department 160 Fieldcrest Avenue Raritan Center Edison, New Jersey 08818	

1
Co

TECHNICAL REPORT DISTRIBUTION LIST, GEN

	<u>No.</u> <u>Copies</u>	
Office of Naval Research Attn: Code 413 800 North Quincy Street Arlington, Virginia 22217	2	Naval Ocean Systems Center Attn: Mr. Joe McCartney San Diego, California 92152
ONR Pasadena Detachment Attn: Dr. R. J. Marcus 1030 East Green Street Pasadena, California 91106	1	Naval Weapons Center Attn: Dr. A. B. Amster, Chemistry Division China Lake, California 93555
Commander, Naval Air Systems Command Attn: Code 310C (H. Rosenwasser) Department of the Navy Washington, D.C. 20360	1	Naval Civil Engineering Laboratory Attn: Dr. R. W. Drisko Port Hueneme, California 93401
Defense Technical Information Center Building 5, Cameron Station Alexandria, Virginia 22314	12	Dean William Tolles Naval Postgraduate School Monterey, California 93940
Dr. Fred Saalfeld Chemistry Division, Code 6100 Naval Research Laboratory Washington, D.C. 20375	1	Scientific Advisor Commandant of the Marine Corps (Code RD-1) Washington, D.C. 20380
U.S. Army Research Office Attn: CRD-AA-IP P. O. Box 12211 Research Triangle Park, N.C. 27709	1	Naval Ship Research and Development Center Attn: Dr. G. Bosmajian, Applied Chemistry Division Annapolis, Maryland 21401
Mr. Vincent Schaper DTNSRDC Code 2803 Annapolis, Maryland 21402	1	Mr. John Boyle Materials Branch Naval Ship Engineering Center Philadelphia, Pennsylvania 19112
Naval Ocean Systems Center Attn: Dr. S. Yamamoto Marine Sciences Division San Diego, California 91232	1	Mr. A. M. Anzalone Administrative Librarian PLASTEC/ARRADCOM Bldg 3401 Dover, New Jersey 07801

END

FILMED

7-83

DTIC